



DESCRIPTION
AXIAL HANDLE, LARGE
T-HANDLE, LARGE
LLIF, TRIAL, 18 X 50 X (8-14mm) -7°
LLIF, TRIAL, 18 X 50 X (8-14mm) -0°
LLIF, PADDLE DISTRACTOR, SIZE: 8-14mm
LLIF, PADDLE SHAVER, SIZE: 8-14mm
LLIF, INSERTER, STRAIGHT
LLIF, BONE GRAFT TAMP
LLIF, GRAFT PACKING BLOCK
LLIF, SLAP HAMMER
MALLET
RETRACTOR BODY
BLADES 90-180mm
SHIM
TROCAR BLADE
SHIM POSITIONER
KEY
TROCAR
DILATOR HOLDER
DILATOR 1&2
SINGLE FIBER OPTIC CABLE
TABLE CLAMP
DISPOSIBLE PEDICLE PROBE
LIGHT MAT

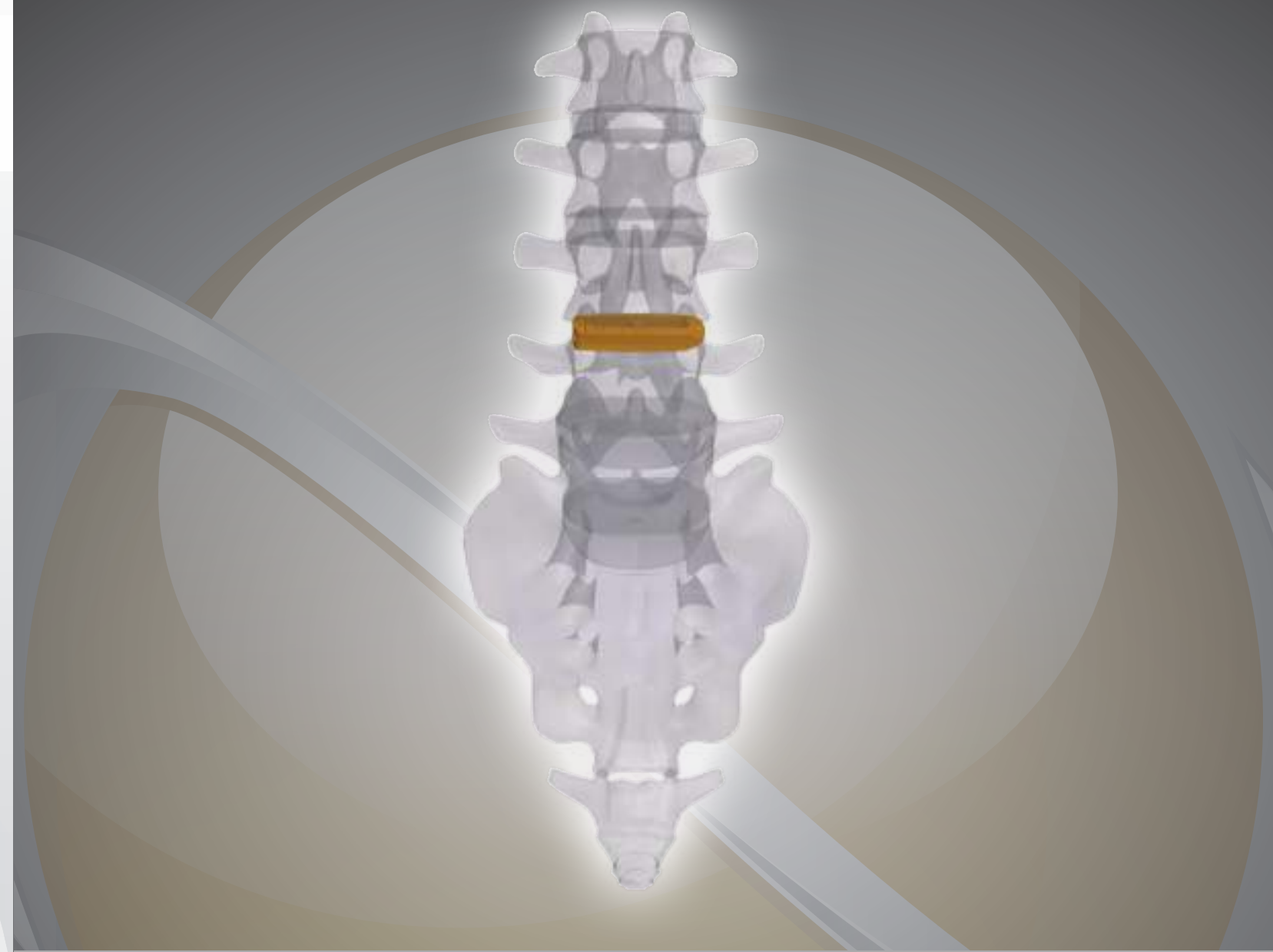

Phone: 760-295-3600
 Fax: 760-295-3610
 Address: 3270 Corporate View, Suite A
 Vista CA 92081
 Email: customerservice@nextorthosurgical.com

EZ-TRACK™
 LATERAL LUMBAR INTERBODY SYSTEM

LLIF

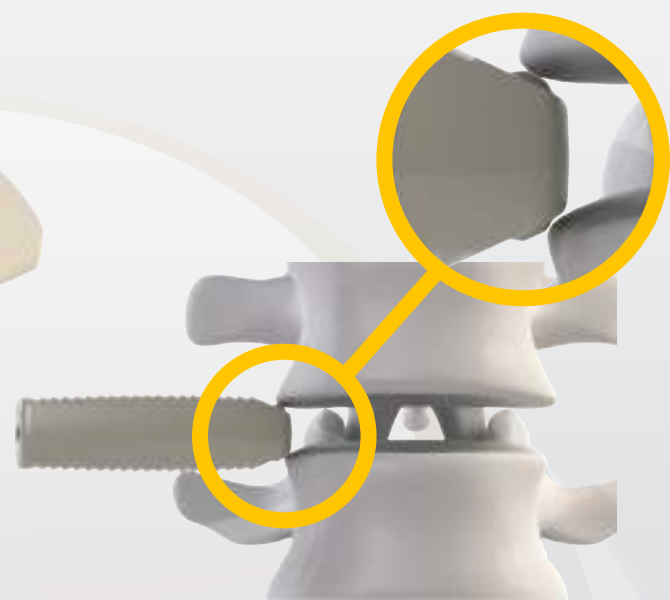

Bone Trap:
The Bone Trap feature retains bone graft throughout insertion.

Tracking Ridge:
Provides precise insertion direction, advanced control and safety.



Biconvex: Contoured to match the vertebral endplates.

LORDOSIS	LENGTH (mm)	WIDTH (mm)	HEIGHT (mm)
0°	40 – 65 (5mm incr.)	18 – 24 (2mm incr.)	8 – 18 (1mm incr.)
7°	40 – 65 (5mm incr.)	18 – 24 (2mm incr.)	8 – 18 (1mm incr.)



Nerve Monitored Dilation:

- 2-step dilation (Ø8mm and Ø14.5mm)
- Nerve monitoring probe



Blades:
Length: 90-180mm



Retractor Body:

- Top and bottom loading blades
- Individual toeing of left and right arms (12.5°)
- Bi-referential table clamp attachments
- Removable handles



Blade Extensions:

- Broach
- Shim